


Data Copy System Instruction for HKPMT precalibration



04 Sep 2026 Kotaro Kanai(UTokyo)
kanai@hep.phys.s.u-tokyo.ac.jp

Schematic Idea

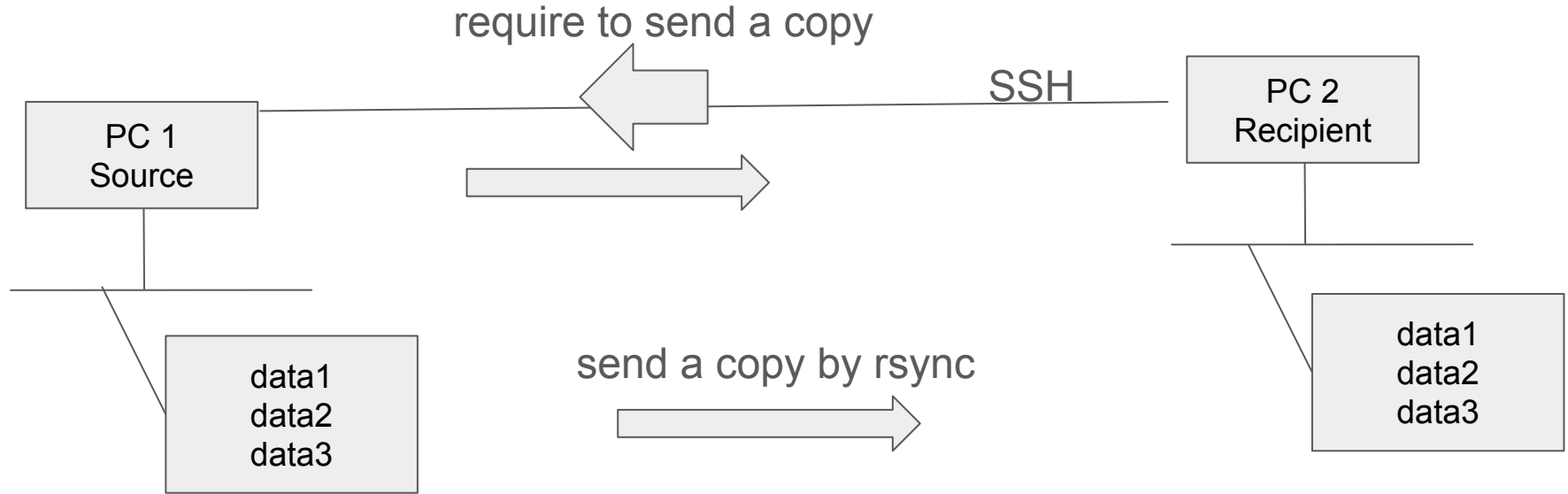

Before execution
You need to create a public key
between PC1 and PC2

You can select

- input folder
- date
- run numbers
- save directory

Before you use it... e.g. PC1=precalib, PC2=your Mac

- You install zenity (GUI tool) to your linux or mac

`sudo apt install zenity`(linux), `brew install zenity`(Mac)

- You have to create a public key in the recipient PC(your Mac)

`ssh-keygen -t ed25519`

↑ key creation. skip 3 questions just push enter key

`ssh-copy-id -p 22 precalib@100.120.120.108`

↑ key sharing with the source PC(precalib) through port 22

Start run

Shell script [0904-cp.sh](https://github.com/kanai0625/precalib-copy) <- download to your working directory in the recipient PC(your Mac)

<https://github.com/kanai0625/precalib-copy>

here suppose you are in the directory [0904-cp.sh](#) exists

Give [0904-cp.sh](#) execution permission ; `chmod +x 0904-cp.sh`

Make your saving directory somewhere you want ; `mkdir your_saving_folder`

Start run in your terminal; `./0904-cp.sh`

input window

PMT ROOT File SCP Transfer Tool

input parameters

Input Folder [initial: data]


Date [initial: 250901]

Run Numbers [e.g., 0 1 2 3 or 0,1,2]

Save Dir [initial: /rawdata/250901]

test files in precalib PC

</home/precalib/Downloads/copytest>

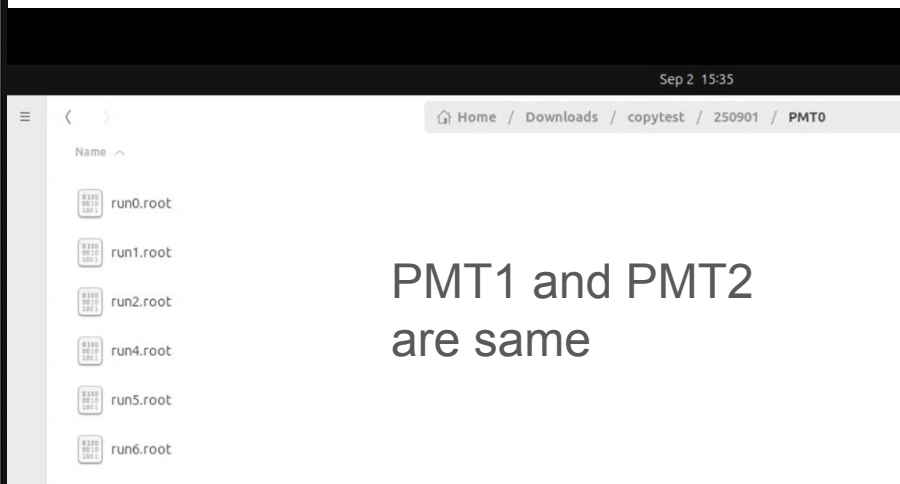

PMT1 and PMT2
are same

Input Folder : data parent directory (basically fixed)


Date : suppose datafiles are in /InputFolder/date

Run numbers : suppose datafile names like run0.root, run1.root etc

Save Dir : your data saving directory

Confirmation window


You confirm your input information and click はい(Y)

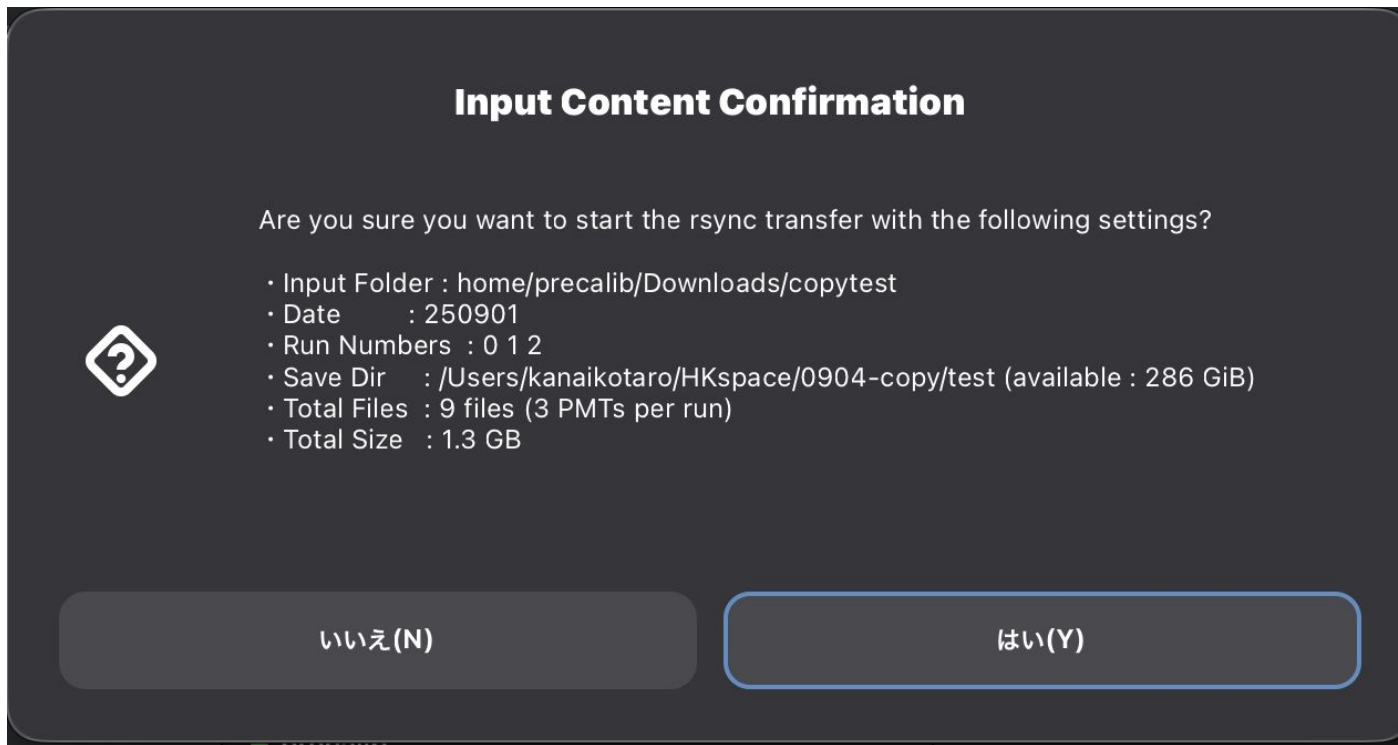

↑ いいえ(No)

↑ はい(Yes)

Progress Window



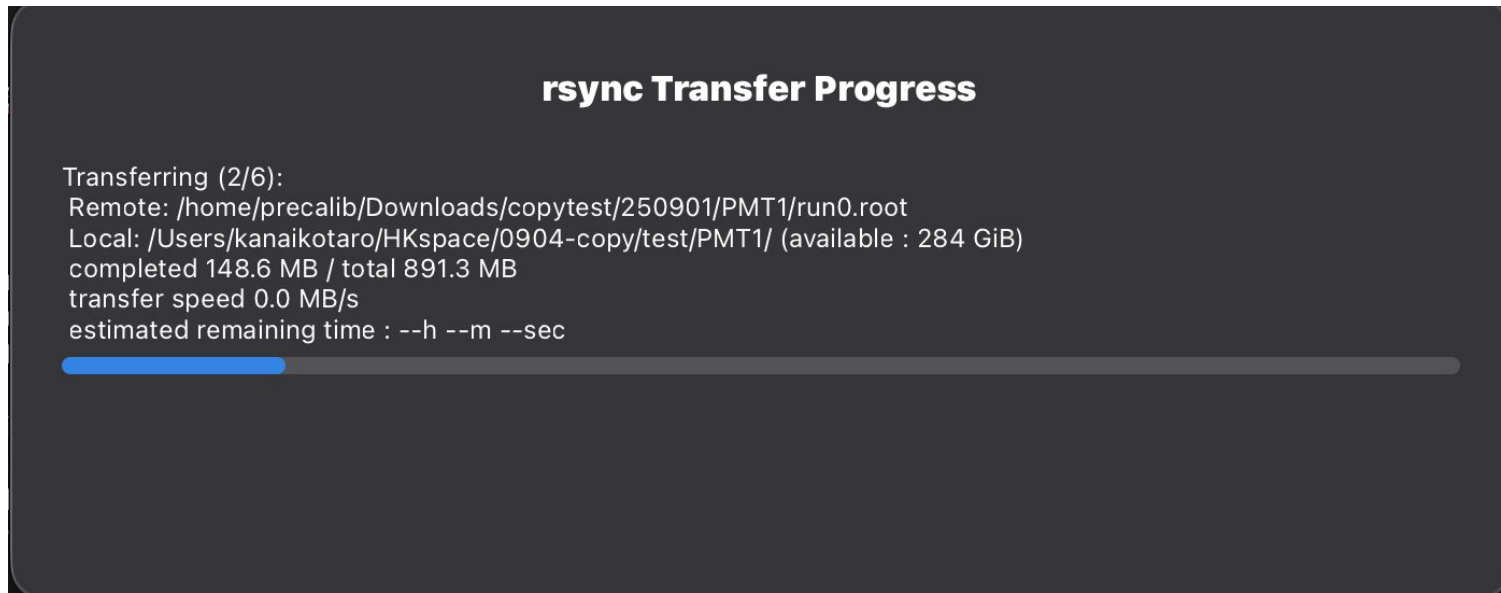


You can see the progress of this copy job realtime

It takes few sec to few min depending on the file size and the number of files u want to copy


Complete Window

When all the files is copied you can see the window

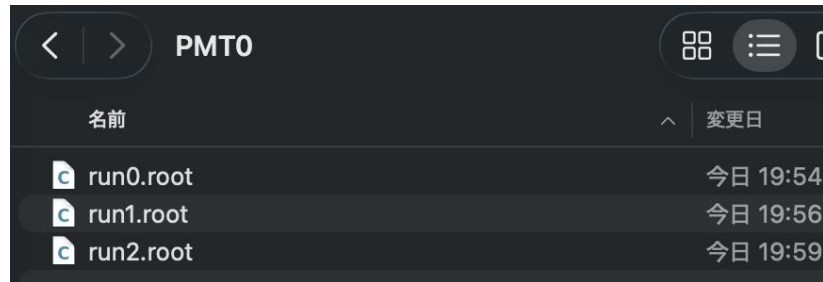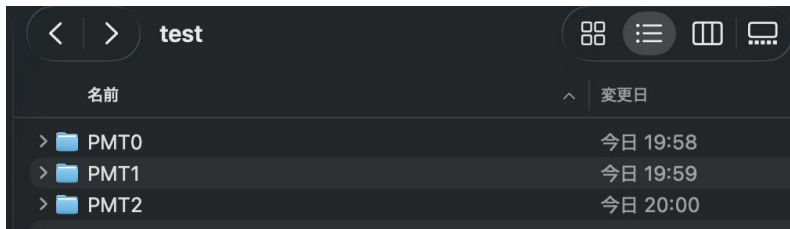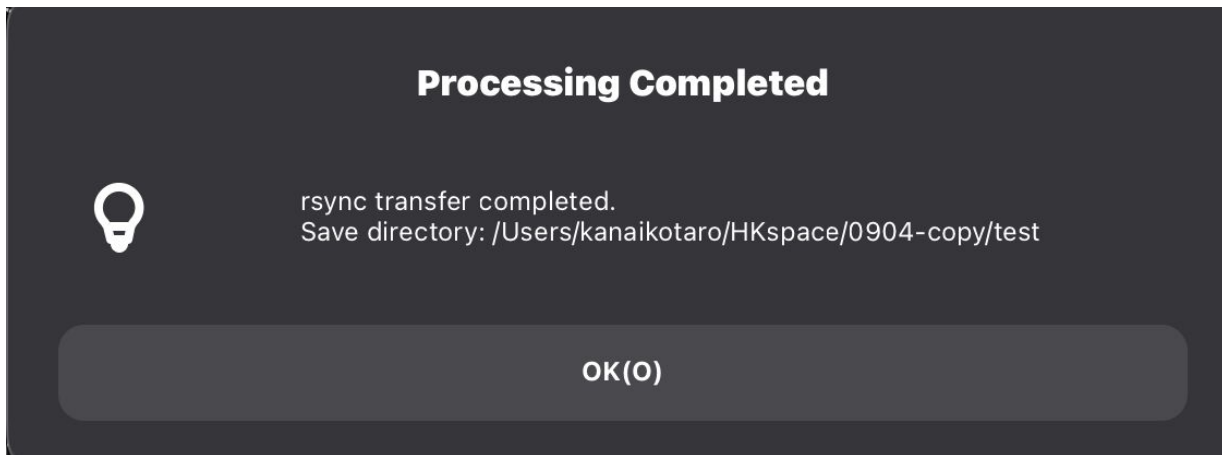

You should confirm that all the files exists in your saving directory

Directory Structure in the source PC 0901version

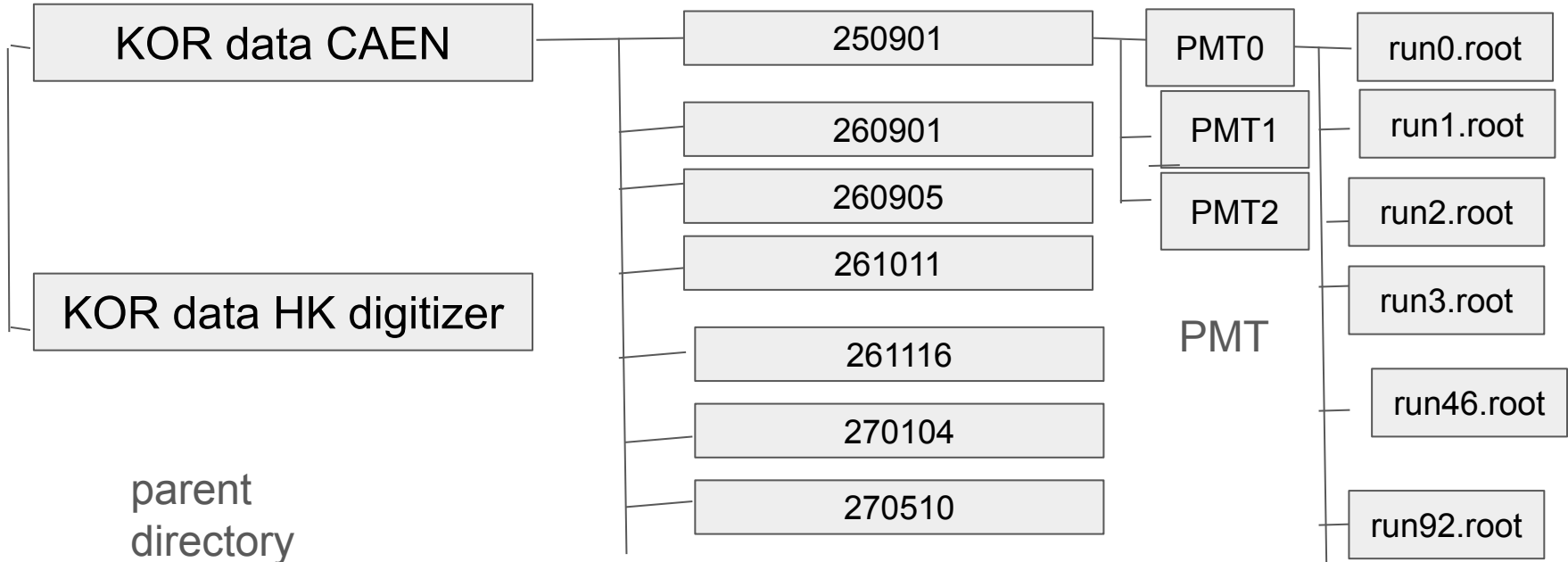

all condition information is supposed to be included in each runN.root file

for later analysis folder does not include experiment variables
You can change the parent directory path modifying 0901-cp.sh

Condition Information Spread Sheet and logbook

Shifters add all conditions they set to measure PMTs like this

date, purpose, time, PMT, HV, B-current, wavelength, digitizer, LASER current etc

Shifters also write these information on a paper logbook kind of a daily report

Data ID	Shift person(s)	Measurement type	Date (JST)	Start Time (JST)	End Time (JST)	PMT (next mon. PMT)	HV [V]	PMT (near door)	HV [V]	B-field current (A)
	-	-	2026/08/24							
	-	-	2026/08/25							
	Junkyo	General scan	2026/08/26							
	Shuoyu	HV scan	2026/08/27	14:43	18:49	EM6400	1630	EL5150	1700	3.12
	Shuoyu	HV scan	2026/08/27	19:04	23:09	EM6400	1730	EL5150	1800	3.12
	Shuoyu	HV scan	2026/08/28	7:17	11:23	EM6400	1830	EL5150	1900	3.12
	Shuoyu	HV scan	2026/08/28	14:18	18:31	EM6400	1880	EL5150	1950	3.12
	Shuoyu	HV scan	2026/08/29	7:22	11:27	EM6400	1580	EL5150	1650	3.12
	Shuoyu	HV scan	2026/08/29	11:40	15:52	EM6400	1480	EL5150	1550	3.12
	Shuoyu	HV scan	2026/08/29	17:44	21:50	EM6400	1480	EL5150	1550	3.12
	Shuoyu	HV scan	2026/08/30			EM6400	1780	EL5150	1850	3.12
	Unik	B-field scan	2026/08/31	18:05	22:15	EM6400	1680	EL5150	1750	3.12
	Unik	B-field scan	2026/08/31	22:30	2:44	EM6400	1680	EL5150	1750	3.15
	Unik	B-field scan	2026/09/01	17:32	21:46	EM6400	1680	EL5150	1750	2.36
	Unik	B-field scan	2026/09/01	21:53	2:07	EM6400	1680	EL5150	1750	1.46
	Kawabata	PMT replacement	2026/09/02							
	Kawabata	Time evolution (short runs)	2026/09/02							
	Kawabata	HV scan	2026/09/03							
	Kawabata	General scan (w/ 1e7 gain)	2026/09/04							
			2026/09/05							
			2026/09/06							
	Yufei	B-field scan	2026/09/07							
	Yufei	B-field scan	2026/09/08							
	Yufei	Temperature scan	2026/09/09							
	Kanai	Temperature scan	2026/09/10							